
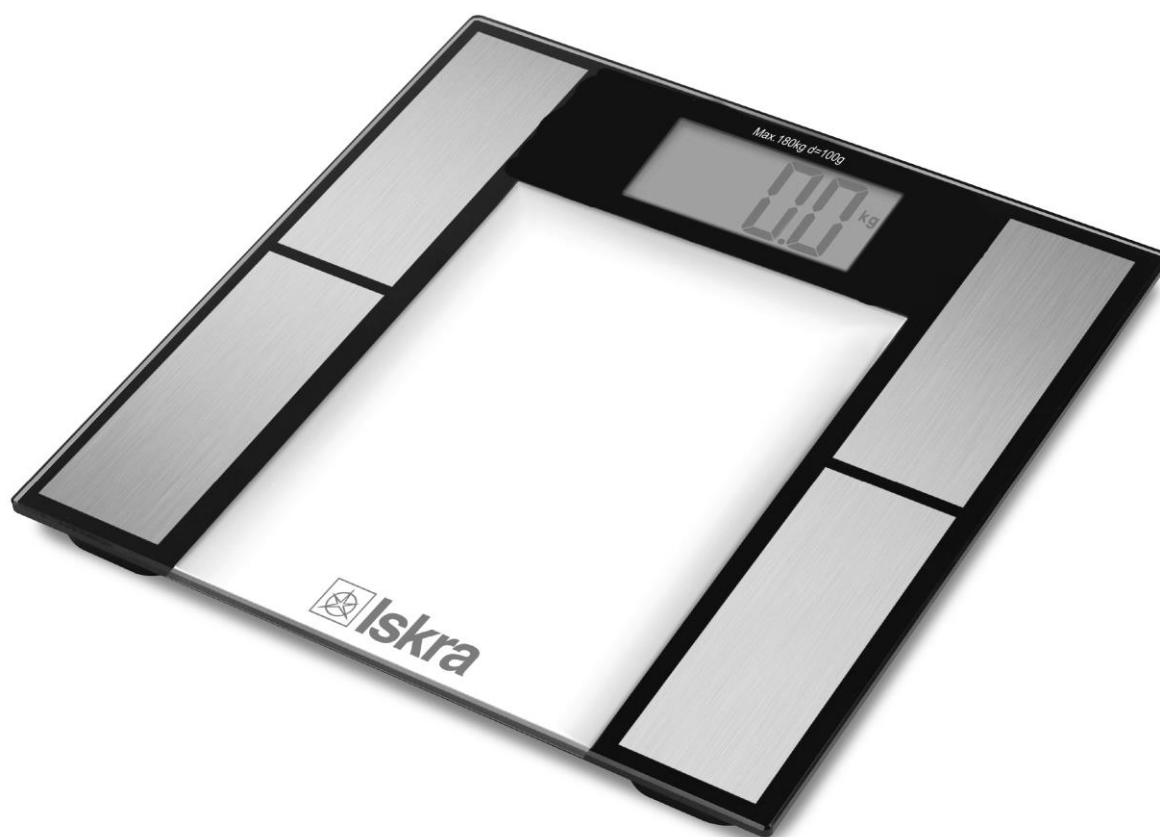
#### GBF1531D-BLE

|         |                                        |                                          |
|---------|----------------------------------------|------------------------------------------|
| SI      | BLUETOOTH® DIAGNOSTIČNA TEHTNICA       | NAVODILA ZA UPORABO                      |
| HR      | BLUETOOTH® DIJAGNOSTIČKA VAGA          | PREVOD ORIGINALNIH UPUTA                 |
| BIH     | BLUETOOTH® DIJAGNOSTIČKA VAGA          | PREVOD ORIGINALNIH UPUTSTVA              |
| SRB/MNE | BLUETOOTH® DIJAGNOSTIČKA VAGA          | PREVOD ORIGINALNIH UPUTSTAVA             |
| MK      | BLUETOOTH® ДИЈАГНОСТИЦИРАЧКА ВАГА      | ПРЕВОД НА ОРИГИНАЛНИ УПАТСТВА            |
| KS      | BLUETOOTH® PESHORE PËR DIJAGNOSTIFIKIM | PËRKTHIMI I UDHËZIMEVE ORIGJINALE        |
| HU      | BLUETOOTH® DIAGNOSZTIKAI MÉRLEG        | AZ EREDETI HASZNÁLATI ÚTMUTATÓ FORDÍTÁSA |
| DE      | BLUETOOTH® DIAGNOSTIKWAAGE             | ÜBERSETZUNG DER ORIGINALANLEITUNG        |
| CZ      | BLUETOOTH® DIAGNOSTICKÁ VÁHA           | PŘEKLAD ORIGINÁLNÍHO NÁVODU              |
| SK      | BLUETOOTH® DIAGNOSTICKÁ VÁHA           | PREKLAD ORIGINÁLNEHO NÁVODU              |
| GB      | BLUETOOTH® DIAGNOSTIC SCALE            | TRANSLATION OF THE ORIGINAL INSTRUCTIONS |

**A**

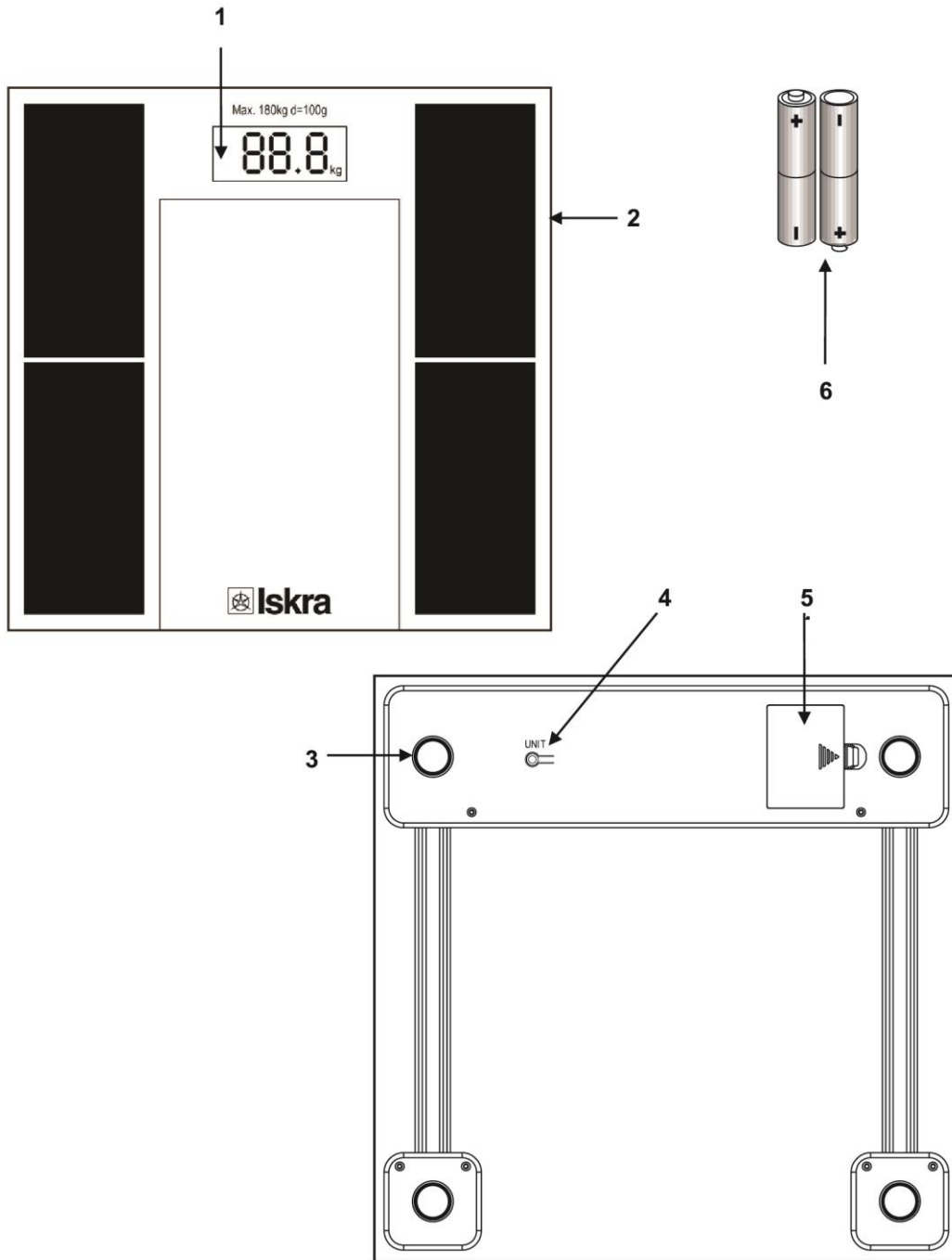**B**

Združljiv z Android4 4.3+  
Bluetooth 4.0+

Kompatibilan s Android4 4.3+  
Bluetooth 4.0+

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**WARNING!**

**BEFORE USING APPLIANCE, CAREFULLY READ THE INSTRUCTIONS WRITTEN IN THE INSTRUCTION MANUAL**

**WARNING!**

**FOLLOW THE INSTRUCTIONS AND SAFETY MEASURES!**

**WARNING**

The manufacturer is not responsible for damage made by improper use, maintenance or modifications of this appliance. Maintenance and modification is limited to the official service personnel.

**WARNING**

Do not throw away this manual. Save it for later use (for example: transfer of ownership)!

**READ AND UNDERSTAND ALL INSTRUCTIONS AND WARNINGS WRITTEN IN THIS MANUAL!**

**CORRESPONDENCE**

In case of mail or telephone correspondence with your salesman or manufacturer, for any reason, concerning your appliance, always specify this information:

- 1) Appliance model-type
- 2) Serial number
- 3) Voltage and frequency
- 4) Date of purchase
- 5) Name of retailer where you bought the appliance
- 6) Precise description of deficiency

**Iskra**

The appliance may only be used by a qualified person who has read this instruction manual and is aware of the dangers connected with using this appliance.

**WARNING**

Failure or refusal to comply with instructions and warnings stated in this manual, can lead to electric shock, fire and/or serious injuries.

In case of using the appliance where metal parts of the appliance may come in contact with hidden electrical lines or its own cable (if installed), always hold the appliance by its plastic isolated handle.

**SAFETY INSTRUCTIONS**

To reduce the chance of fire, electric shock or personal injury, follow the instructions written in this manual, when using this appliance.

## **USAGE IN ACCORDANCE WITH REGULATIONS**

The manufacturer is not responsible for the damage caused by improper use of the tool or by unauthorised modifications made by the user.

Follow all safety, assembly and usage instructions.

## **IMPORTANT INSTRUCTIONS**

This appliance is not intended for commercial use. Save these instructions for future reference. Before you start using your appliance, it's very important to read and follow the instructions, which will also come in handy in the future.

Pay special attention to the chapter which deals with your personal safety.

### **For your safety**

- This appliance is not intended for use by persons (including children under the age of 8) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should not be allowed to play with packing material. Risk of asphyxiation.
- Do not overload the appliance.
- Check the appliance before each use. If the appliance is damaged in any way, do not use the appliance.
- Do not submerge the appliance in water or other liquid.
- Only use this appliance indoors. Do not expose this appliance to rain.
- Store this appliance so it doesn't fall into the bath tub or sink.
- Do not insert any objects inside the appliance.
- Place the appliance on a leveled, dry and stable surface.
- Do not expose the appliance to high temperatures (e.g. stove) or weather conditions (e.g. rain). Do not pour liquid into the appliance.
- Do not use the appliance if its housing is damaged, broken or deformed.
- Do not change or modify the appliance in any way. Only use the appliance the way it is intended.
- Be careful not to drop the appliance. Do not hit the appliance or expose it to strong vibrations or magnetic fields.

- When storing this appliance for a longer period of time, do not put any weight on its surface.

## BATTERY HANDLING

### **WARNING!**

**If children or animals swallow the battery, they can seriously endanger their health.**

- Store batteries out of reach of children and animals.
- Discard used batteries and batteries leaking battery acid. Discard batteries in accordance with regulations.
- If children or animals swallow the battery, seek medical help immediately.

## PRODUCT SPECIFIC WARNINGS

### **WARNING!**

**Do not use this appliance if you are dependent on medical implants (if you have for example a cardiac pacemaker). Measuring voltage for bioelectrical impedance analysis may affect the operation of such implants.**

- Do not use this appliance during pregnancy.
- Do not step on the appliance with wet feet. Make sure the appliance is dry, before you step on it.
- Remove footwear, before stepping on the appliance.
- Do not step on the edge of the appliance.
- Never jump on the appliance or let it fall on the ground, when moving it.
- Clean the appliance with a wet cloth. Do not submerge the appliance into water.
- Do not expose the appliance to heat or extreme cold.

## PARTS DESCRIPTION (Figure A)

- |                         |                                           |
|-------------------------|-------------------------------------------|
| 1. LCD display          | 4. Switching between units of measurement |
| 2. Measuring electrodes | 5. Battery compartment                    |
| 3. Rubberized legs      | 6. Batteries 2x LR03/R03/AAA              |

**Supplied with two 1.5V LR03/R03/AAA batteries.**

## LCD DISPLAY (Figure B)

- |                              |                               |                       |
|------------------------------|-------------------------------|-----------------------|
| 1. Age                       | 8. Total body water (TBW) (%) | 15. Normal            |
| 2. Bone mass %               | 9. Percentage display         | 16. Underfat          |
| 3. Numerical display         | 10. Body fat percentage       | 17. Personal data     |
| 4. Body height (foot/inch)   | 11. Body mass index           | 18. Female            |
| 5. Body height (centimeter)  | 12. Basal metabolic rate      | 19. Male              |
| 6. Body weight (kilogram)    | 13. Obese                     | 20. Muscle percentage |
| 7. Body weight (stone/pound) | 14. Overfat                   |                       |

## MODE OF OPERATION

This scale uses a harmless and undetectable low level electric current for measuring. This current runs through the body, where it encounters different resistance. The electrical resistance is then measured and calculated together with other body properties (height, age and sex) to obtain the percentage of water, bone and muscle. This manner of measuring is called »Bioelectrical impedance analysis« (BIA).

In case of certain body properties or deficiencies, measured data may deviate from normalised values. This applies to persons under the age of 10 years and to persons with significantly short or long legs. Data may deviate with persons with fever, osteoporosis, illnesses that cause swelling, patients with dialysis and persons that take cardiac medications.

Measurements performed by this appliance are merely approximation from a medical point of view. Exact measurements of body fat, water, muscle and bone mass can only be performed by a medical specialist with specialised medicinal equipment.

## USAGE

### Before first use

- Remove all packaging
- Remove the battery cover on the bottom side of the appliance and remove the plastic stripe between the battery and contact.
- Close the battery compartment, so the cover locks audibly.

## SWITCHING BETWEEN MEASURING SYSTEMS

By pressing the button on the bottom of the appliance, you can switch between different measuring units: kg / lb / st:lb.

## MEASURING WEIGHT

### WARNING!

**Do not overload the appliance. Maximum weight of the appliance is: 180 kg.**

### WARNING!

**Do not step on the appliance with wet feet. Make sure the appliance is dry before stepping on it.**

- Place the appliance on levelled surface.
- Step on the appliance carefully. Stand in the centre of the appliance and keep straight posture.
- Your measured weight will show on the display.

For the most accurate results, weigh every day at approximately the same hour.

## MEASURING BODY FAT, WATER, MUSCLE AND BONE MASS

- Stand barefoot on the appliance and keep straight posture. You must keep your feet on the appliance the whole time, to maintain contact with electrodes.
- Your weight will display first. The appliance will then analyse the measurements provided by the entered data and measured weight.
- The measurements will display in a rotating order: body weight (**kg**), body fat percentage (**fat%**), total body water (**TBW%**), muscle percentage (**mus%**), bone mass (**bone kg**), body mass index (**BMI**) and basal metabolic rate (**KCAL**).

### SCALE USE WITHOUT APPLICATION

If you use the scale without application, the only correct measurement is your weight. All other measurement values will be calculated using standard parameters in the case of use without application, therefore the readings can be more or less inaccurate.

The measured data is not stored in the scale, if you use it without an application.

**We recommend using the scale with the application.**

### SCALE AND APPLICATION (Figure C)

The scale is intended for weighing the body and analyzing body data. The measured data is sent via a Bluetooth® connection to a smartphone or tablet. With free OKOK International application (hereinafter referred to as "application") you can store this data on a smartphone or tablet and analyze it. Data is always visible only to the current user, but not to other users.

### COMPATIBILITY

Scale and application are compatible with the following smartphones and tablets with Bluetooth® Smart Ready (from version 4.0 onwards):

- Apple® iPhone® 4s and higher
- Apple® iPad® 3rd generation and later versions
- Apple® iPad mini™ 1st generation and later versions
- Apple® iPad Air™ and later versions
- Smartphones and tablets with Android™ 4.3 or higher

A more detailed list of compatible devices can be found on the following website:

[http://www.tookok.cn/download\\_en.html](http://www.tookok.cn/download_en.html)

### CREATING USER ACCOUNT

The scale has unlimited user saving slots in which you and your family members save your settings and measurement values. For this purpose, it is necessary for each user to create their own user account.

In addition to personal information you will also need to enter: gender, date of birth and height. These data will be used as parameters for calculating relative values, e.g. for BMI.

1. To specify a user in OKOK International, follow the step-by-step instructions to register.
2. Select the symbol of the "Me" application to perform user-specific or product-specific settings.

### ESTABLISH A Bluetooth® CONNECTION

1. Make sure you have Bluetooth® enabled on your phone or tablet.
2. Run the application and the "Home" menu will open. 
3. Select the "Device" icon in the application.
4. To establish a connection between the scale and the smartphone, or tablet select the "Device" and press "Pair" to enter scale nickname.
5. The scale will be switched on.  
Each user should step on the scale once with bare feet, so that it is completely determined as a user of the scale.

### ADDING FAMILY MEMBERS IN APPLICATION

To set up more family members for the scale, follow these steps:

1. Open the application.
2. Select the "Me" symbol in the application.
3. Then select "Family Member".
4. Then press "Add".
5. Enter the user name and the new user data.
6. After starting the initial measurements, setting up a new user is completed.

### NOTE

**When weighing family members, note that you manually assign measurements to the appropriate user from the list in the application.**

## EVALUATION OF MEASUREMENT RESULTS

### NOTE

**Measuring will be precise only if you input correct height, age and gender data.**

Changes in measured values depend heavily on sport activity and diet. When you analyse measured results it is most sensible to observe long term changes in body water, fat and muscle values. Higher muscle values and lower weight and fat values are signs of well-balanced diet or training.

### NOTE

**The tables bellow are for reference only. Consult your doctor for exact information.**

#### Body mass index (BMI) - WOMEN

| Age   | Underweight | Normal | Overweight | Obese |
|-------|-------------|--------|------------|-------|
| 16-24 | < 19        | 19-24  | 25-28      | > 28  |
| 25-34 | < 20        | 20-25  | 26-29      | > 29  |
| 35-44 | < 21        | 21-26  | 27-30      | > 30  |
| 45-54 | < 22        | 22-27  | 28-31      | > 31  |
| 55-64 | < 23        | 23-28  | 29-32      | > 32  |
| > 65  | < 24        | 24-29  | 30-33      | > 33  |

#### Body mass index (BMI) - MEN

| Age   | Underweight | Normal | Overweight | Obese |
|-------|-------------|--------|------------|-------|
| 16    | < 19        | 19-24  | 25-28      | > 28  |
| 17-24 | < 20        | 20-25  | 26-29      | > 29  |
| 25-34 | < 21        | 21-26  | 27-30      | > 30  |
| 35-44 | < 22        | 22-27  | 28-31      | > 31  |
| 45-54 | < 23        | 23-28  | 29-32      | > 32  |
| 55-64 | < 24        | 24-29  | 30-33      | > 33  |
| > 65  | < 25        | 25-30  | 31-34      | > 34  |

#### Body fat percentage WOMEN

| Age   | Underfat | Normal | Overfat | Obese |
|-------|----------|--------|---------|-------|
| 10-29 | < 23     | 23-34  | 35-45   | > 45  |
| 30-59 | < 23     | 23-34  | 35-45   | > 45  |
| 60-69 | < 23     | 23-34  | 35-45   | > 45  |
| > 70  | < 23     | 23-37  | 38-45   | > 45  |

#### Body fat percentage MEN

| Age   | Undefat | Normal | Overfat | Obese |
|-------|---------|--------|---------|-------|
| 10-29 | < 11    | 11-24  | 25-35   | > 35  |
| 30-59 | < 13    | 13-25  | 26-35   | > 35  |
| 60-69 | < 12    | 12-25  | 26-35   | > 35  |
| > 70  | < 12    | 12-29  | 30-35   | > 35  |

#### Total body water WOMEN

| Age   | Low  | Normal | High |
|-------|------|--------|------|
| 10-29 | < 52 | 52-72  | > 72 |
| 30-49 | < 40 | 40-60  | > 60 |
| > 50  | < 45 | 45-65  | > 65 |

### Total body water MEN

| Age   | Low  | Normal | High |
|-------|------|--------|------|
| 10-29 | < 55 | 55-75  | > 75 |
| 30-49 | < 50 | 50-70  | > 70 |
| > 50  | < 48 | 48-68  | > 68 |

### Muscle percentage WOMEN

| Age   | Low  | Normal | High |
|-------|------|--------|------|
| 10-29 | < 35 | 35-45  | > 45 |
| 30-59 | < 28 | 28-40  | > 40 |
| 60-69 | < 23 | 23-35  | > 35 |
| > 70  | < 23 | 23-30  | > 30 |

### Muscle percentage MEN

| Age   | Low  | Normal | High |
|-------|------|--------|------|
| 10-29 | < 40 | 40-50  | > 50 |
| 30-59 | < 33 | 33-45  | > 45 |
| 60-69 | < 32 | 32-40  | > 40 |
| > 70  | < 30 | 30-35  | > 35 |

### SHUT-OFF

- The appliance will shut off automatically after a few seconds of being idle.

### WARNING INDICATORS

- If the display shows »Lo« indicates that the battery is empty. In this case, replace the battery on the bottom of the appliance.
- If the display shows »Err«, indicates that the appliance is over encumbered.
- If the display shows »Err2«, indicates a fat measuring error.
- If the measured body fat percentage is lower than 5%, the display will show »FATL« indication.
- If the measured body fat percentage is higher than 50%, the display will show »FATH« indication.

### CLEANING AND MAINTENANCE

### WARNING!

**Do not submerge the appliance into water or other liquid.**

- Do not use abrasive cleaning accessories to clean this appliance.
- Clean the appliance with a moist cloth and let it dry completely.
- This appliance is not intended for commercial use.
- Store this appliance in a dry space out of reach of children.

### SERVICE

Contact an authorized service agent for any kind of repairs.

A damaged power cable should be replaced with original parts by an authorized service agent.

**Before handing the appliance over to the service agent, make sure it's cleaned and in its original packaging.**

| Model:       | Maximum weight: | Power supply:                             |
|--------------|-----------------|-------------------------------------------|
| GBF1531D-BLE | 180 kgs         | 2 x 1.5 V DC battery<br>type LR03/R03/AAA |

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